



TRIBUTE™
MANUFACTURED BY SOLARIS

Therapeutic NightWear™ Garments

Distributed in Canada by
**LYMPHEDEMA
DEPOT LTD**
www.lymphedemadepot.com

Tribute Therapeutic NightWear Garments

Long-Term Self-Management

Tribute custom designed NightWear garments are essential for a successful lymphedema home management program. Worn nightly, your Tribute garment will offer an effective solution to help regulate fluid fluctuations, decrease tissue turgor, and assist a weakened lymphatic system. Solaris' attention to the details provides you with a garment that eliminates complex bandaging, and is stress free, comfortable, easy to care for, and backed with a one year warranty.



Table of Contents

Tribute Components & Warranties	4
Helpful Hints	5
Washing Instructions	6
Tribute Arm Garments	8
Tribute Leg Garments	10
Tribute Torso Garments	12
Tribute Head Garments	14
Outer Jackets, VCJs & Contraindications	15
Contact Information	16

Tribute Components

- Meryl / Lycra fibers
- Latex-free
- Medical grade foam
- Manufactured in the USA

Your Warranty

Tribute garments, made with our proprietary fabric, are warranted to be free of material defects for one year. Solaris will repair or replace any garment found to be defective. Using bleach, fabric softener or dryer sheets, as well as any other additives, will void the warranty. Exposure to heat over 130° F will damage the garment and void the warranty. Damage due to fingernail tears, animal damage, urine saturation, or excessive pulling are not covered by the warranty. Solaris does not warrant Outer Jackets or Variable Compression Jackets.

Helpful Hints

- Inside out is right! Seams are on the exterior of the garment to prevent skin irritation.
- Take your time when putting on your Tribute. Use care and be gentle to prevent thread breakage.
- To prevent nail tears, wear rubber gloves when donning and doffing.
- During the day, store your Tribute in the mesh bag provided.
- Wait at least 10 minutes to put on your Tribute after applying skin moisturizer.
- For more Helpful Hints, watch the Tribute product videos on the Solaris website: www.solarismed.com.

Washing Instructions

Tribute Garments

Machine wash your Tribute in warm water on the permanent press cycle, using a regular detergent. To extend the life of your garment, double rinse to ensure all soap residue is removed. We also suggest washing your Tribute with a normal load of laundry.

Dry on permanent press. Heat exceeding 130° F will permanently damage the fabric, deteriorate the foam filling, and void your warranty.

Do not use fabric softener in either the washing or drying cycle. Using any form of additives, such as bleach or fabric softener, will void your warranty.

We suggest washing your garment prior to donning it for the first time.

Washing Instructions cont.

Outer Jackets & Variable Compression Jackets

Hand wash Outer Jackets (OJs) & Variable Compression Jackets (VCJs) as frequently as needed to eliminate dirt and odor.

Dip in warm water with a mild detergent, rinse thoroughly and gently squeeze out excess water.

Lay flat to dry. Machine washing and drying may shorten the useful life of this garment.

OJs & VCJs should be replaced every 6 months, depending on your wearing protocol, and are not covered under the Tribute warranty.


Tribute Arm Garments

Easy On:

Step 1 Fold the upper half of the garment over the lower half.

Step 2 Gently draw the garment over your hand (the sewn seam should align with the thumb).

Step 3 Roll the upper half of the arm sleeve onto the arm, twisting the garment slightly toward your body to align the sewn seam to the center of the armpit.

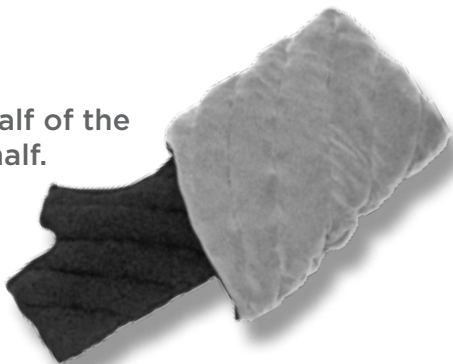

Tribute Arm Garments

Easy Off:

Step 1 Gently unroll the upper half of the garment onto the lower half.

Step 2 Remove the garment completely.

Step 3 Store during the day in the mesh bag provided to allow airflow and the evaporation of any moisture.

(Do not wash your Tribute in the mesh bag)



Tribute Leg Garments

Easy On:

Step 1 Fold the upper half of the garment over the lower half.

Step 2 Carefully pull the garment over your foot using the heel assist loop.

Step 3 Adjust for comfort:

- a. Chevron - the sewn seams are positioned along the inside and outside of the leg.
- b. Vertical - the sewn seams are positioned along the front & back of the leg.

Step 4 Roll the upper half of the garment onto your leg. Do not twist the garment as you unroll it.



Tribute Leg Garments

Easy Off:

Step 1 Gently unroll the upper half of the garment onto the lower half.

Step 2 Remove the garment completely.

Step 3 Store during the day in the mesh bag provided to allow airflow and the evaporation of any moisture.

(Do not wash your Tribute in the mesh bag)



Tribute Torso Garments

Easy On:

Step 1 Position the garment on the body, with the zipper in the front.

Step 2 Zip the garment slowly and carefully.

Step 3 When applicable, adjust shoulder straps as needed to enhance the fit while standing or sleeping.

Easy Off:

Unzip the garment slowly and carefully, store in mesh bag.

(Do not wash your Tribute in the mesh bag)



Tribute Head Garments

Easy On:

Step 1 Position the unit so that the *shorter strap* will wrap around the neck.

Step 2 Using buckles and Velcro, attach neck straps lightly.

Step 3 Attach longer straps at the crown of the head. Adjust the straps to provide a comfortable fit.

(For more detailed instructions, please see your Mandibular Unit user guide.)



Outer Jackets & VCJs

OJs - Easy On/Off:

Gather your Outer Jacket the length of the garment. Place the OJ over your hand or foot, then gently slide it up the Tribute garment.

VCJs - Easy On/Off:

Your therapist will adjust the straps to fit your limb. Be sure not to separate the VCJ from the Tribute when donning or doffing. Simply zip and unzip.


Contraindications

- Arterial Insufficiency or Degeneration
- Deep Vein Thrombosis
- Untreated Congestive Heart Failure or Cancer
- Untreated Localized or Systemic Infection
- Absent or Severely Impaired Sensation

Adjusting Finger Spacers

Adjusting Finger Spacers

To adjust the length of the finger spacers, turn the garment inside out. Simply snip the finger spacers a little bit at a time to decrease the length.

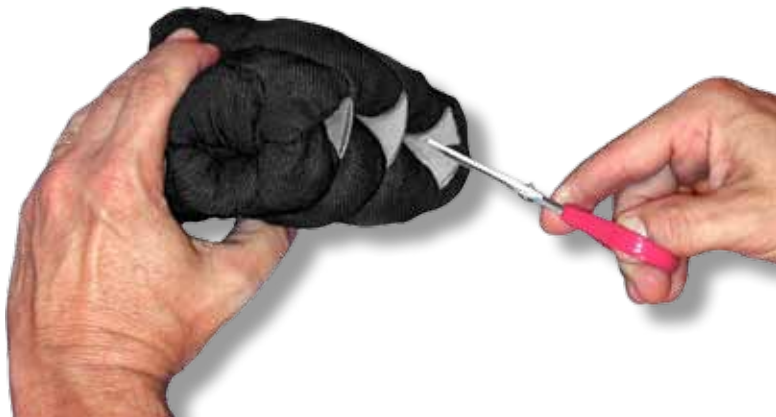

SOLARIS
FREEDOM TO LIVE



Contact Lymphedema Depot
T: (905) 687-8500
F: (905) 687-8224
E: info@lymphedemadepot.com

www.lymphedemadepot.com



SOLARIS 6737 West Washington Street Suite 3260 West Allis WI 53214

custom
made

11/01/10